


RBS-330

USER'S MANUAL



330lb x 0.2lb /
150kg x 0.1kg

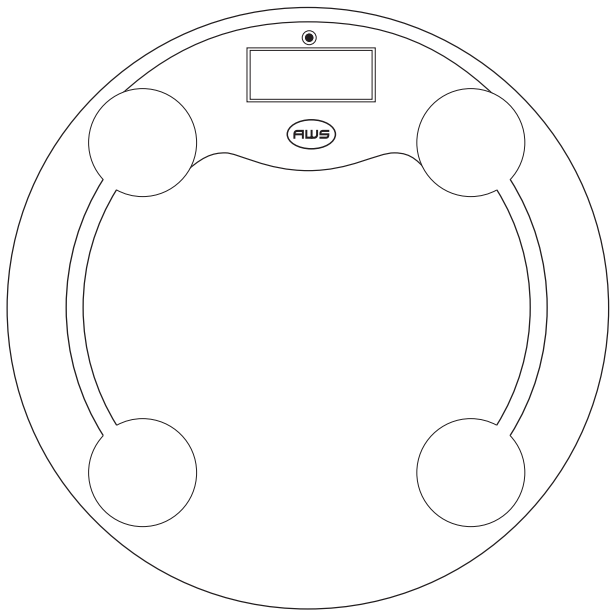


AMERICAN WEIGH
SCALES INC.

American Weigh
RBS-330 User's Manual
330lb x 0.2lb

Thank you for purchasing the RBS-330 Round Bathroom Scale. With normal care and proper treatment, it will provide you with years of reliable service. Please read all operating instructions carefully before using the scale.

- Avoid any exposure to extreme heat or cold. Your scale will work better when operating in normal room temperature. Always allow the scale to acclimate to normal room temperature for at least an hour before use.
- Give the scale sufficient warm up time before use. Usually 30 to 60 seconds to give the internal components time to stabilize.
- Keep your scale in a clean environment. Dust, dirt, moisture, vibration, air currents, or a close proximity to other electronic equipment can cause an adverse effect on the reliability and accuracy of your scale.
- Handle with care. Gently apply all items to be weighed onto the weighing surface. Although this scale is designed to be quite durable, avoid rough treatment as this may permanently damage the internal sensor.
- Do not shake, drop, or otherwise shock the scale. This is a precision instrument and must be handled with extreme care.
- Only operate the scale on a stable, vibration free surface.
- Always store in a dry place to avoid moisture accumulation on internal components.

Parts Description

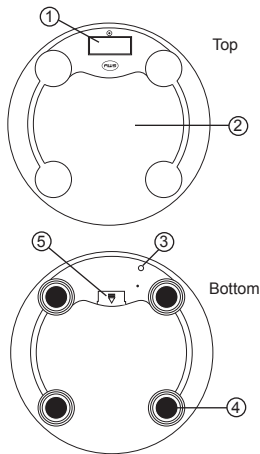

1. LCD Display
2. Tempered Glass
3. LB/KG Switch
4. Scale foot
5. 9V Battery Compartment

I. Setup

Remove the main unit from the package. Insert a 9V battery (included) into the battery compartment (see above). You are now ready to weigh!

II. Weighing

Scale may not turn on if weight is applied too slow. In this case, you can tap the scale at any time to turn the display on.

1. The scale is ready to weigh at any time.
3. Step onto the scale and wait for the readout to stabilize.
4. The weight will be displayed on the LCD
5. After 10 seconds, the scale will turn off (Auto-Off)
6. Replace the battery if “LO” is displayed on the LCD.

III. Mode Selection (See prev. page for button locations)

1. Flip the scale upside down and locate the LB/KG button.
2. Press the LB/KG button to switch between metric or imperial weighing modes. (Note the change on the LCD)
3. Set the scale right-side-up again when mode is selected.

IV. Specifications

Capacity:	330lb ; 150kg
Graduation:	0.2lb; 0.1kg
Units:	pounds (lb) ; kilograms (kg)
Display:	4 digits 1” in height